

OraGRAFT® Endure

Moldable Demineralized Fibers with Cancellous



The OraGraft Advantage

✓ OPTIMIZED HANDLING

OraGraft Endure was developed specifically for grafting procedures where there is a need for a bone graft material that has osteoinductive potential, maintains space and can easily be formed to fit a defect. OraGraft Endure is a mix of the demineralized moldable fibers with cancellous particulate in a 1:1 ratio. This mix maximizes the amount of cancellous particulate and the amount of moldable demineralized fibers without compromising handling.


Feature	Benefit
Hydrates quickly	Saves chair time
Moldable	Conforms and stays contained in the surgical site



✓ HOSPITABLE SCAFFOLD

The cortical fibers in OraGraft Endure create a hospitable scaffold, and the cancellous particulate provides the added benefit of space maintenance.

Feature	Benefit
Long cortical fibers with multiple protrusions and a rough surface	Provides a large surface area and many contact points for cellular attachment ¹
Interconnected cortical fibers	Allows cells to easily spread out and make connections with each other ¹
Favorable porosity and pore size ²⁻¹²	Promotes cellular attachment and proliferation
Lot tested with nude mouse assay	Confirmation that final product has osteoinductive potential
Cancellous particulate 250-1000 microns	Provides increased space maintenance compared to fibers alone

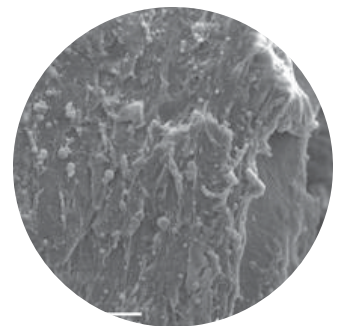

✓ 100% BONE, NO CARRIER


OraGraft Endure is comprised of 100% human bone. The cortical fibers have been demineralized using PAD® technology to encourage natural remodeling during the bone-healing process.

Feature	Benefit
Demineralized using PAD technology	Optimal osteoinductive potential ¹³⁻¹⁷
No carrier	100 percent bone ^{1,18}
Human bone	Facilitates natural remodeling during the bone healing process

Offering a 1:1 ratio of cancellous particulates with moldable demineralized fibers, OraGraft Endure is designed for enhanced osteoinductive potential with a scaffold shown to maintain space and support angiogenesis. The combination has been shown to be effective in supporting guided bone regeneration with the primary application use for treating fenestration defects.¹⁹⁻²¹

Demineralized using proprietary PAD technology that targets optimal residual calcium levels of 1-4% without compromising the graft's inherent osteoinductive potential.

Processed using Allowash XG technology, which provides a Sterility Assurance Level (SAL) of 10^{-6} while protecting the graft's inherent biomechanical and biochemical function.

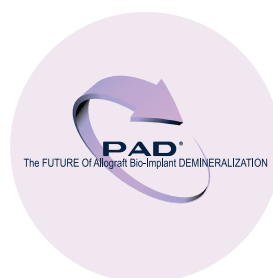

REFERENCES

- Murphy MB, Suzuki RK, Sand TT, et al. Short term culture of mesenchymal stem cells with commercial osteoconductive carriers provides unique insights into biocompatibility. J Clin. Med. 2013; 2:49-66; doi:10.3390/jcm2030049
- Data on file at LifeNet Health, DHF 14-006 (Porosity Profile of Compressed Demineralized Bone Matrix Fiber Strips Protocol ES-15-076, 11/17/2015)
- Karageorgiou, V. and D. Kaplan, Porosity of 3D biomaterial scaffolds and osteogenesis. Biomaterials, 2005. 26(27): p. 5474-91.
- Kuboki, Y., et al., BMP-induced osteogenesis on the surface of hydroxyapatite with geometrically feasible and nonfeasible structures: topology of osteogenesis. J Biomed Mater Res, 1998. 39(2): p. 190-9.
- Eggl, P.S., W. Muller, and R.K. Schenk, Porous hydroxyapatite and tricalcium phosphate cylinders with two different pore size ranges implanted in the cancellous bone of rabbits. A comparative histomorphometric and histologic study of bony ingrowth and implant substitution. Clin Orthop Relat Res, 1988(232): p. 127-38.
- Flautre, B., et al., Porous HA ceramic for bone replacement: role of the pores and interconnections - experimental study in the rabbit. J Mater Sci Mater Med, 2001. 12(8): p. 679-82.
- Galois L, et al., Bone ingrowth into two porous ceramics with different pore sizes: An experimental study. Acta Orthop Belg, 2004. 70(6): p. 598-603.
- Gauthier, O., et al., Macroporous biphasic calcium phosphate ceramics versus injectable bone substitute: a comparative study 3 and 8 weeks after implantation in rabbit bone. J Mater Sci Mater Med, 2001. 12(5): p. 385-90.
- Hulbert, S.F., et al., Potential of ceramic materials as permanently implantable skeletal prostheses. J Biomed Mater Res, 1970. 4(3): p. 433-56.
- Motomiya, M., et al., Effect of Hydroxyapatite porous characteristics on healing outcomes in rabbit posterolateral spinal fusion model. Eur Spine J, 2007. 16(12): p. 2215-24.
- Shimazaki, K. and V. Mooney, Comparative study of porous hydroxyapatite and tricalcium phosphate as bone substitute. J Orthop Res, 1985. 3(3): p. 301-10.
- Uchida, A., et al., Bone ingrowth into three different porous ceramics implanted into the tibia of rats and rabbits. J Orthop Res, 1985. 3(1): p. 65-77.
- Zhang M, Powers RM, and Wolfenbarger L. Effect(s) of the demineralization process on the osteoinductivity of demineralized bone matrix. J Periodontol. 1997; 68:1085-1092
- Turonis JW, McPherson JC 3rd, Cuenin MF, et al. The effect of residual calcium in decalcified freeze-dried bone allograft in a critical-sized defect in the Rattus norvegicus calvarium. J Oral Implantol. 2006;32(2):55-62
- Herold RW, Pashley DH, Cuenin MF, et al. The effects of varying degrees of Allograft Decalcification on the Cultured Porcine Osteoclast cells. J Periodontol. 2002 Feb; 73(2):213-9
- Mott DA., Mailhot J, Cuenin MF, Sharawy M, Borke J. Enhancement of osteoblast proliferation in vitro by selective enrichment of demineralized freeze-dried bone allograft with specific growth factors. J Oral Implantol 2002;28(2):57-66
- Pietrzak WS, Ali SN, Chitturi D, Jacob M, Woodell-May JE. BMP depletion occurs during prolonged acid demineralization of bone: characterization and implications for graft preparation. Cell Tiss. Bank 2007 (Published online)
- Lee JH, Baek H-R, Lee KM, et al. The effect of poloxamer 407-based hydrogel on the osteoinductivity of demineralized bone matrix. Clinics in Orthopedic Surgery. 2014; 6(4):455-461; doi:10.3390/jcm2030049
- Rosen, PS, Reynolds, MA. Guided Bone Regeneration for Dehiscence and Fenestration Defects on Implants Using an Absorbable Polymer Barrier, J. Periodontol. Feb. 2001 Vol. 72, No. 2, 250-256
- Fugazzotto PA, Shanaman R, Manos T, Shectman R, Guided Bone Regeneration Around Titanium Implants: Report of the Treatment of 1,503 Sites With Clinical Reentries, The International Journal of Periodontics & Restorative Dentistry, Vol. 17, No. 3 1997, 293-299
- Eskow AJ, Mealey BL. Evaluation of Healing Following Tooth Extraction With Ridge Preservation Using Cortical Versus Cancellous Freeze-Dried Bone Allograft, J. Periodontol April 2014, Vol. 85, No. 4. 514-524



ariston dental
better answers